ELICOM
ELECTRONIC

VAAKATALO

Vaakatalo Oy
Vestonkatu 11
33580 TAMPERE
info@vaakatalo.com
www.vaakatalo.com

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Si-10

Si-10S

Si-10SP

CONTENTS

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RIGHTS AND LIABILITIES

All rights reserved.

No patent liability is assumed with respect to the use of the information contained herein. While every precaution has been taken in the preparation of this book, Elicom Electronic assumes no responsibility for errors or commissions. Neither is any liability assumed for damages resulting from the use of the information contained herein.

The information herein is believed to be both accurate and reliable. Elicom Electronic however, would be obliged to be informed if any errors occur. Elicom Electronic cannot accept any liability for direct or indirect damages resulting from the use of this manual.

Elicom Electronic reserves the right to revise this manual and alter its content without notification at any time.

Neither Elicom Electronic nor its affiliates shall be liable to the purchaser of this product or third parties for damages, losses, costs, or expenses incurred by purchaser or third parties as a result of: accident, misuse, or abuse of this product or unauthorized modifications, repairs, or alterations to this product, or failure to strictly comply with Elicom Electronic cooperating and maintenance instructions. Elicom Electronic shall not be liable against any damages or problems arising from the use of any options or any consumable products other than those designated as Original Elicom Electronic products.

NOTICE: The contents of this manual are subject to change without notice.

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3.2 PERIPHERAL CONNECTIONS (ENG)

Warning:

1. *Please always remember that Si30 indicator is very low voltage measuring instrument. Please do not forget that the instrument must be powered off before inserting or removing any peripheral connector.*
2. *The quality of instrument's grounding will provide weighing accuracy and the safety of your indicator. If the energy condition of your plant is bad, prepare a special power line and grounding.*
3. *All required electrical connections should be done as described below.*

3.2.1 MAIN BOARD WIRING CONNECTIONS

To avoid damages, the load cell wiring should be made carefully before energizing the instrument. Load cell connection detail is shown in Kuva 3..

In 4-wire installations jumpers RJ4 and RJ5 must be shortened. The two jumpers are located right from the load cell connector.

If you have junction box in your system, use 6 wire cable between indicator and the junction box.

In 6-wire installation you must remove RJ4 and RJ5

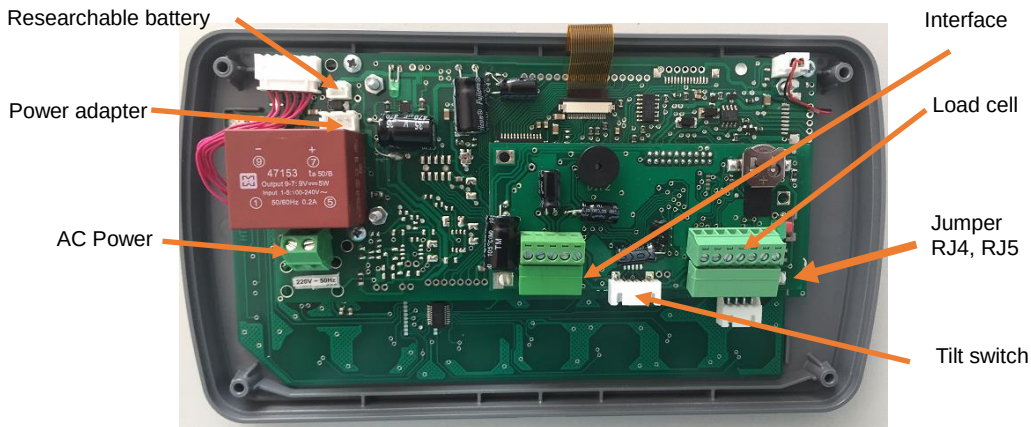
Figure 3.2.1.1 Main board connections AC model

The electrical connectors of I10 are located at the bottom of the indicator as seen below left and electrical terminals of I10S are located on the printed circuit board in the housing as seen below right.

3.2.2 EXTERNAL CONNECTIONS

To avoid damages, the load cell wiring should be made carefully before energizing the instrument. Load cell connection, Interface and power detail is shown in Kuva 3.

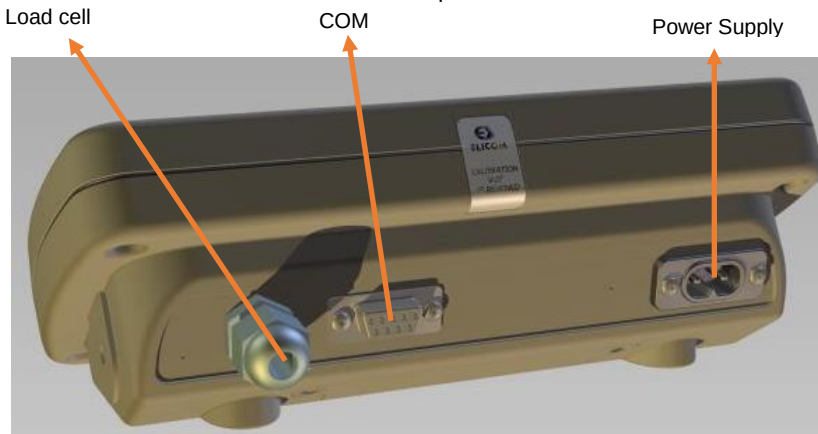
Figure 3.1 External connection Si10 AC model



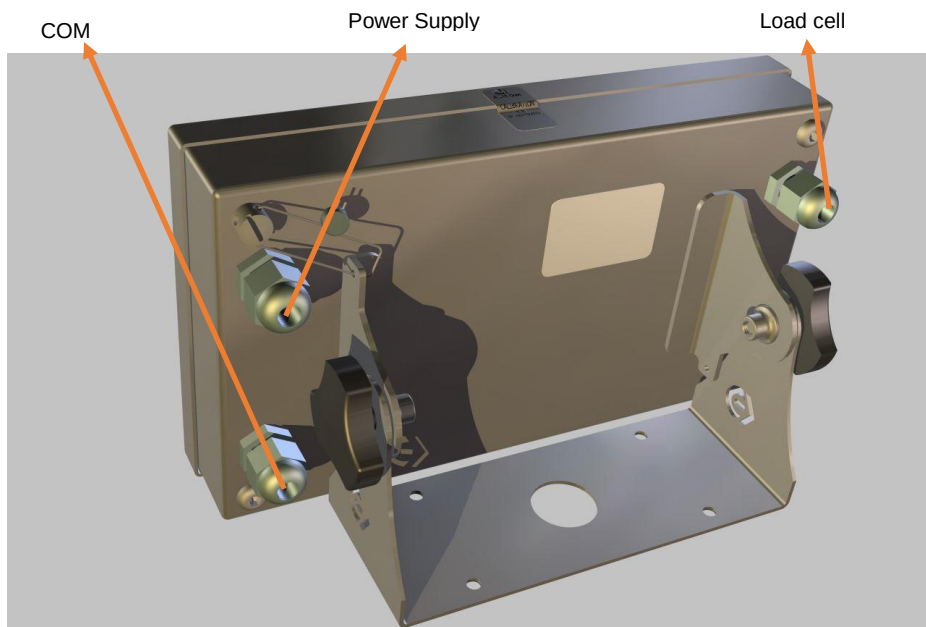
Figure 3.2 External connection Si10S AC model

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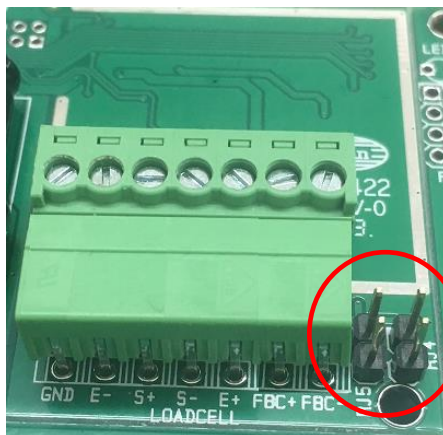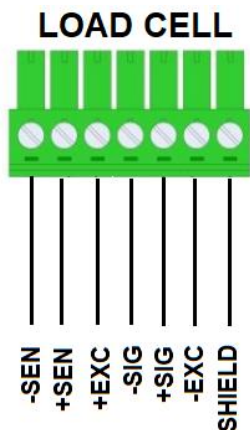
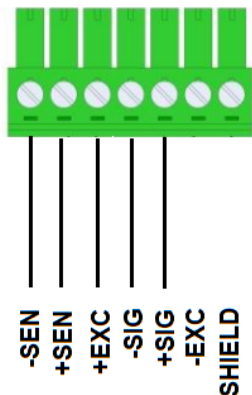
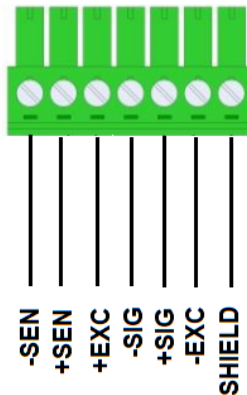
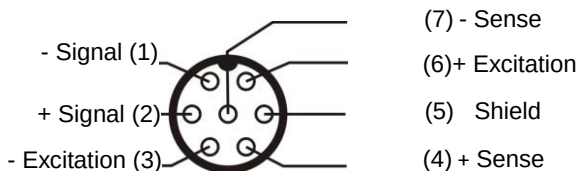
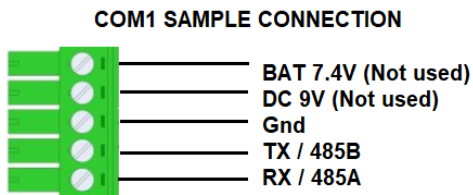
Figure 3.5 Sample COM1 internal termination connection

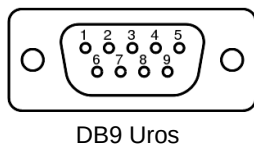	Male DB9	
Signal	COM1	COM2
TXD	2	8
RXD	3	7
GND	5	5

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
4.13 TILT SWITCH /EI SAATAVILLA

The instrument can be equipped with an optional external tilt switch.

When the instrument is out of the Tilt limits the switch will be active and weighing will be terminated by the following messages:


TILTED SCALE	Tilt Switch is ON
TILT SENSOR	Problem with Tilt Sensor
TILT SETUP	Tilt sensor not found or problem in calibration data

5. ASETUSVALIKKO KAAVIO

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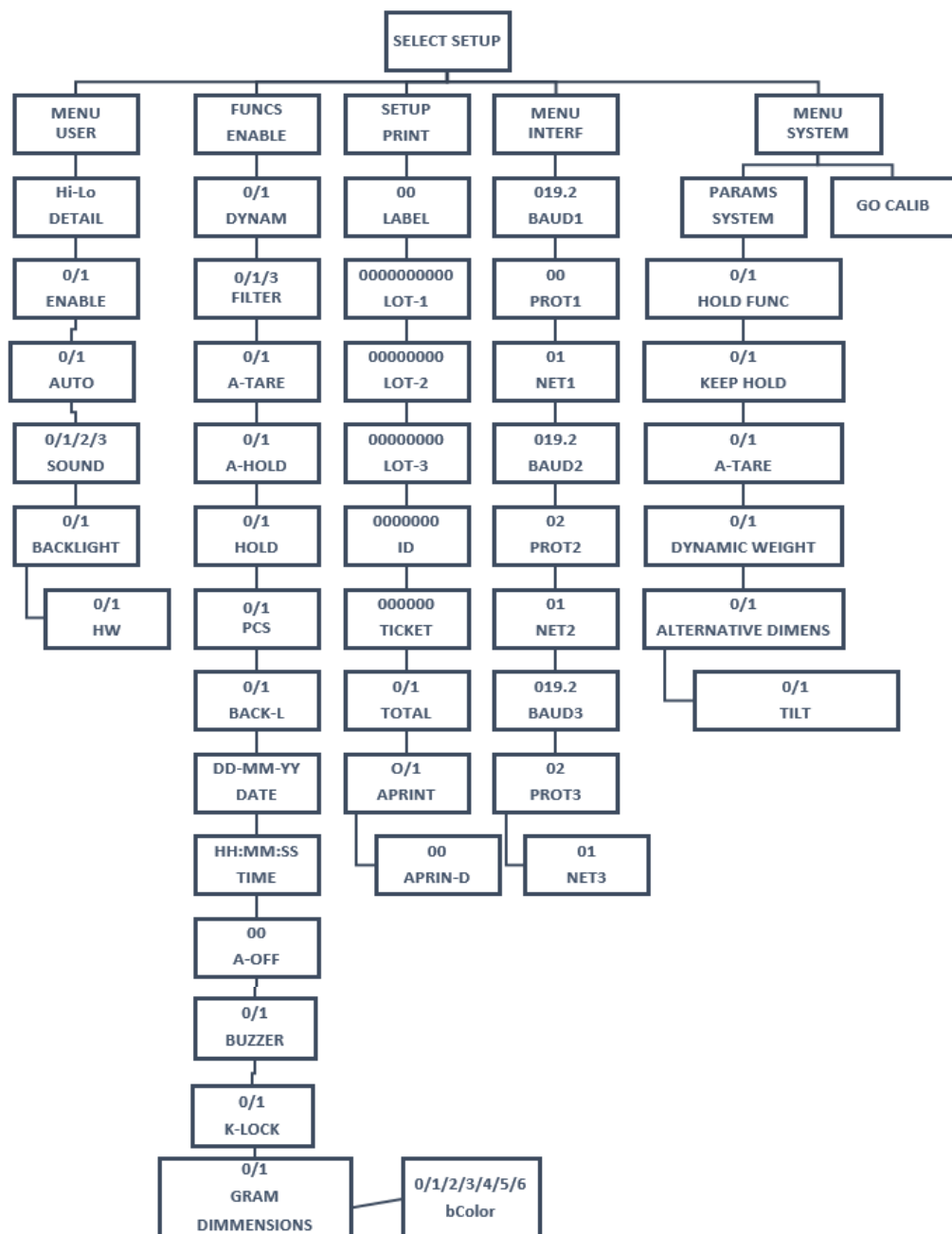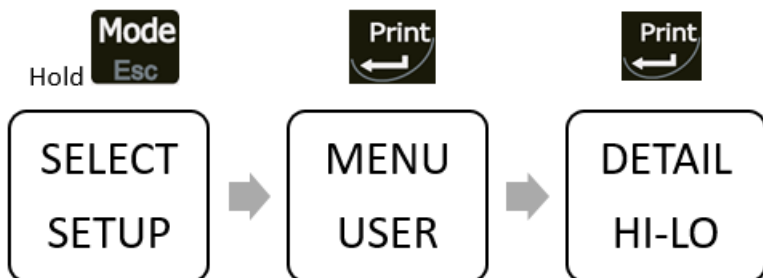6.1 CHECKWEIGHING FUNCTION ACTIVATION

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To change press

Press to go to next parameter

To exit menu hold until indicator is in weighing mode

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
10 INTERFACE SETUP (ENG)

Setup menu for setting up interfaces
Access from general menu as shown:

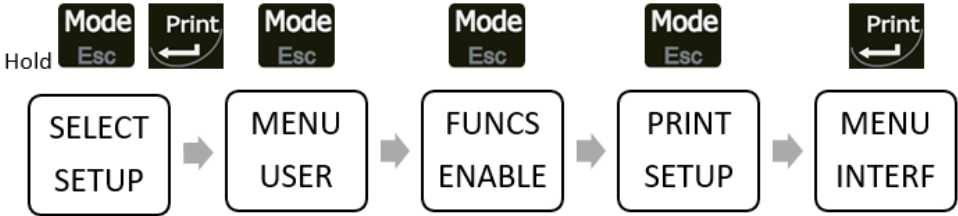10.1 COM1 SPEED

Enter COM1 speed:

004.8 – 4800bps

009.6 – 9600bps

019.2 - 19200bps (used for PC software)

Use arrows to input

Press  to go to next parameter

10.2 COM1 PROTOCOL

Select protocol

00- PC software SiComm

01- Label editor connection

02- Printer connection

03- Reserved

04- Reserved

05- Reserved

Use arrows to input

Press  to go to next parameter

10.3 COM1 SUBPROTOCOL

Used for protocol variations

00- Reserved

01- Reserved

02- Reserved

03- Reserved

04- Reserved

Use arrows to input

Press  to go to next parameter

019.2

BAUD1

00

PROT-1

0

NET-1

10.4 COM2 SPEED

Enter COM1 speed:

004.8 – 4800bps

009.6 – 9600bps

019.2 - 19200bps (used for PC software)

Use arrows to input

Press to go to next parameter

10.5 COM2 PROTOCOL

Select protocol

00- PC software SiComm

01- Label editor connection

02- Printer connection

03- Reserved

04- Reserved

05- Reserved

Use arrows to input

Press to go to next parameter

10.6 COM2 SUBPROTOCOL

Used for protocol variations

00- Reserved

01- Reserved

02- Reserved

03- Reserved

04- Reserved

Use arrows to input

Press to go to next parameter

10.7 COM3 SPEED

(Reserved for internal devices)

Enter COM3 speed:

004.8 – 4800bps

009.6 – 9600bps

019.2 - 19200bps

Use arrows to input

Press to go to next parameter

**019.2
BAUD2**

**00
PROT2**

**00
NET-2**

**019.2
BAUND3**

10.8 COM3 PROTOCOL

Select protocol

- 00- PC software SiComm
- 01- Label editor connection
- 02- Printer connection
- 03- Reserved
- 04- Reserved

Use arrows to input

Press  to go to next parameter

10.9 COM3 SUBPROTOCOL

Used for protocol variations

- 00- Reserved
- 01- Reserved
- 02- Reserved
- 03- Reserved
- 04- Reserved

Use arrows to input

Press  to go to next parameter

To exit menu hold  until indicator is in weighing mode

IMPORTANT: To save interface setting you must go

over all parameters with  until the instrument makes a click signal and exits Interface Menu

00
PROT-3


00
NET-3

11 SYSTEM PARAMETERS (LUKITTU, ENG)

System parameters menu is jumper protected and it is used during initial scale setup.

From this menu you can enable/disable different functions

Access from general menu as shown:

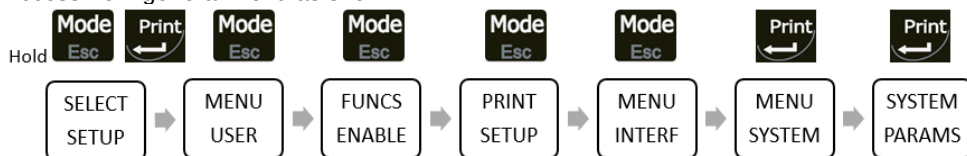11.1 FUNCTIONS HOLD AND AUTO HOLD

- 0- Disabled
- 1- Active

To change press

Press to go to next parameter

11.2 CLEAR HOLD/AUTO HOLD

- 0- Clears hold result after 8secs
- 1- Keep hold result until manually cleared with C key

To change press

Press to go to next parameter

Note: In legal for trade scales this parameter must be set to "0"

11.3 AUTO TARE

- 0- Disabled
- 1- Active

To change press

Press to go to next parameter

0
HOLD FUNCS

0
KEEP HOLD

0
A-TARE FUNC

11.4 DYNAMIC WEIGHING

- 0- Disabled
- 1- Active

To change press

Press to go to next parameter

11.5 UNIT CHANGE

Allows unit change

- 0- Disabled
- 1- Active

To change press

Press to go to next parameter

11.6 TILT SWITCH

The instrument can be equipped by internal or external TILT switch depending on application.

- 00- Disabled
- 01- Disabled
- 02- Disabled
- 03- Active internal tilt switch (if equipped)
- 04- Active external tilt switch active

Use numerical keyboard to input

Press to go to next parameter

If Tilt Switch is active the following menus appear

0
DYNAMIC WEIGHT

0
ALTERN DIMENS

0
TILT TYPE

11.7 FLAT POSITION CALIBRATION

This parameter is used to calibrate the TILIT switch when scale in flat position (leveled).
Top line is showing the current weight.
Level the scale and place load.

To set the level calibration point hold 
Message GET POSIT will appear if calibration is made successfully

Press  to go to next parameter

11.8 TILT POSITION CALIBRATION

This parameter is used to calibration the TILIT switch when scale in tilted position (not leveled).
Top line is showing the current weight.
Set the scale to position where the on the TITL switch should be active scale and place load.

To set the level calibration point hold 
Message GET POSIT will appear if calibration is made successfully

Press  to go to next parameter

To exit menu hold  until indicator is in weighing mode

0.000
FLAT POSIT

0.000
TILTED POSIT

Main office

3 Kapitan Krystev str.

Silistra 7500

BULGARIA

Tel.: +359 86 820 200, 882 255

e-mail: info@elicom-bg.com

www.elicom.bg

Vaakatalo Oy
Vestonkatu 11
33580 TAMPERE
info@vaakatalo.com
www.vaakatalo.com

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Copyright © 2017 by Elicom Electronic